



www.adipogen.com


Progranulin – Marker of Neuroinflammation

Progranulin (PGRN) is a highly conserved secreted cysteine-rich protein that is expressed in multiple cell types, both in the CNS and in peripheral tissues. PGRN is composed of seven ~6kDa granulin (GRN) proteins and acts both directly and via its conversion to granulins, to regulate multifunctional biological activities, including major roles in cancer, inflammation, metabolic disease and neurodegeneration. PGRN is an abundant, non-conventional, stress-induced, extracellular matrix-bound secreted growth factor-like molecule and chaperone. PGRN binds to several functionally different receptor families in a cell/tissue-specific and condition/disease-dependent manner. For example, PGRN binding with TNFR and DR3 has an important anti-inflammatory role in immune cells, particularly Tregs and macrophages. PGRN/Ephrin type-A receptor 2 (EphA2) interaction is involved in the proliferative influence of PGRN. PGRN binds and activates Notch receptors, enhancing the regenerative capacity of injured neurons. PGRN is also a lysosomal resident protein. Sortilin and lipoprotein receptor-related protein 1 (LRP1) have been demonstrated to be the lysosomal trafficking receptors for PGRN with the help of Prosaposin. In the brain, PGRN plays a key role in the development, survival, function, and maintenance of neurons and microglia. It regulates lysosomal biogenesis, inflammation, repair, stress response and aging. Absence of progranulin in microglia causes increased production and release of multiple cytokines, suggesting that PGRN regulates microglia activation. PGRN seems to affect microglial proliferation, recruitment, differentiation, activation and phagocytosis, suggesting that PGRN plays a central role in the regulation of neuroinflammatory responses. In neurons, PGRN i) co-localizes in late endosomes and early lysosomes with the transmembrane protein TMEM106B, ii) co-localizes with markers such as BDNF along axons, iii) influences synaptic structure and function at synaptic and extrasynaptic sites, where it is secreted in an activity-dependent manner, and iv) extracellular PGRN is endocytosed and delivered to lysosomes. The lysosomal function of PGRN is not well characterized, but probably involves regulation of proteins such as cathepsins, glucocerebrosidase (GBA) or TMEM106B and likely contributes to neurodegener-

ation (see Figure). PGRN loss-of-function mutations cause neuronal ceroid lipofuscinosis or frontotemporal Lobar Dementia (FTLD) in a gene dosage-dependent manner and PGRN is therefore **a valuable biomarker of FTLD**. Mutations that reduce PGRN levels also increase the risk for developing Alzheimer's disease, Parkinson's disease and other neurological diseases.

SELECTED REVIEWS: The lysosomal function of progranulin, a guardian against neurodegeneration: D.H. Paushter, et al.; Acta Neuropathol. 136, 1 (2018) • Progranulin: A conductor of receptors orchestra, a chaperone of lysosomal enzymes and a therapeutic target for multiple diseases: Y. Cui, et al.; Cytokine Growth Factor Rev. 45, 53 (2019) • Progranulin as a therapeutic target in neurodegenerative diseases: H. Rhinn, et al.; Trends Pharmacol. Sci. 43, 641 (2022) • Lysosomal functions of progranulin and implications for treatment of frontotemporal dementia: M.J. Simon, et al.; Trends Cell Biol. (Epub ahead) (2022)

FIGURE: The roles and binding partners of progranulin in neuronal cells.

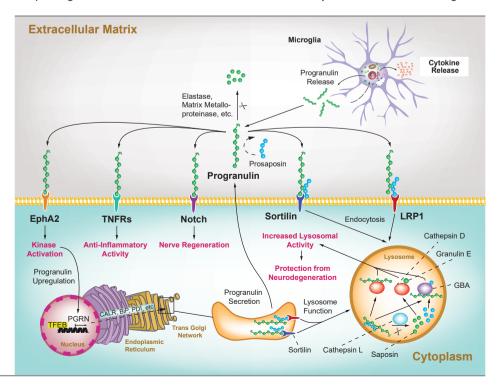

NEW & Unique Progranulin ELISA Kit Assays

Measuring human GRN mutations in Frontotemporal Diseases

GRN mutations are frequent causes of familial frontotemporal lobar degeneration (FTLD). Progranulin expression in plasma predicts GRN mutations status, independently of symptom onset proximity. Progranulin loss-of-function mutations are among the most frequent genetic causes, responsible for 20% of familial FTLD. Progranulin levels in plasma or serum, constitute a reliable, cost-effective biomarker, suitable as a screening tool in patients with familial frontotemporal degeneration.

AdipoGen's Progranulin (human) ELISA Kits are the only and unique kits on the market detecting all mutations of Progranulin. They are both validated on Mutated and Null *GRN* plasma samples.





NEW mAb/mAb-based human Progranulin ELISA Kit

Progranulin (human) ELISA Kit (mAb/mAb)

AG-45B-0027 96 wells

Sensitivity: 60 pg/ml

Range: 0.063 ng/ml – 4 ng/ml

Sample: Cell Culture Supernatant, Plasma, Serum

Specificity: Detects human Progranulin.


The new mAb/mAb based human Progranulin ELISA Kit has been thoroughly validated and compared to the standard pAb/pAb based ELISA Kit (Prod. No. AG-45A-0018Y). With this new ELISA Kit, levels below 50 ng/ml are strongly suggestive of *GRN* mutations. In a validation on 191 patient samples, confirmed by a molecular gene analysis, the new kit provided a sensitivity and specificity of 100% for detecting FTLD mutations.


Concentrations of Progranulin plasma levels in FTLD and healthy samples:

Using the new mAb-based ELISA kit (AG-45B-0027), plasma progranulin levels in FTLD-GRN are lower (~28 ng/ml) compared to controls (~98 ng/ml).

ROC curve of the new ELISA Kit:

ROC analysis shows that at a threshold of 55 ng/ml plasma progranulin measured with ELISA Kit (AG-45B-0027) discriminates FTLD-GRN patients from controls with 100% sensitivity and 100% specificity.

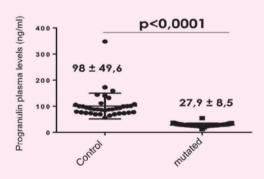

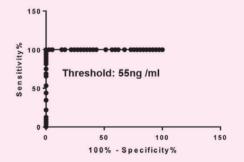

The Standard pAb/pAb-based human Progranulin ELISA Kit

Progranulin (human) ELISA Kit

AG-45A-0018YEK 96 wells AG-45A-0018YTP 2 x 96 wells

Sensitivity: 32 pg/ml

Range: 0.063 ng/ml – 4 ng/ml

Sample: Cell Culture Supernatant, Plasma, Serum, Urine

Specificity: Detects human Progranulin.


LIT: Plasma progranulin levels for frontotemporal dementia in clinical practice: a 10-year French experience: L. Sellami, et al.; Neurobiol. Aging 91, 167.e1 (2020)

The Standard pAb/pAb-based human Progranulin ELISA Kit has been widely used by the research community and thoroughly validated on Mutated and Null *GRN* plasma samples.

- Trusted Reproducible Results!
- Used to Determine Cut-Off values for FTLD!
- Cited in Hundreds of Scientific Publications!

Values measured by the AG-45A-0018Y (pAb-based) & AG-45B-0027 (mAb-based) ELISA Kits:

A strong positive correlation was observed between the pAb-based ELISA kit (Reference kit, AG-45A-0018Y) and the newly developed mAb-based ELISA kit (AG-45B-0027) in the measurement of Progranulin plasma levels.

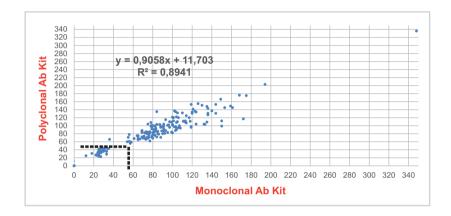

The Standard Mouse & Rat Progranulin ELISA Kits

| ASSAYS | PID | SIZE | SAMPLE |
|-------------------------------|--------------|-------------------------|---------------------------------|
| Progranulin (mouse) ELISA Kit | AG-45A-0019Y | 96 wells 2 x 96 wells | Cell Culture Supernatant, Serum |
| Progranulin (rat) ELISA Kit | AG-45A-0043Y | 96 wells 2 x 96 wells | Cell Culture Supernatant, Serum |

Progranulin Antibodies & Tagged Proteins

| ANTIBODIES | PID | SIZE | ISOTYPE/SOURCE | APPLICATION | SPECIES |
|-------------------------------|--------------|-----------------|----------------|----------------|---------|
| anti-Progranulin (human), pAb | AG-25A-0112 | 100 µg | Guinea pig | ELISA, IHC, WB | Hu |
| anti-Progranulin (mouse), pAb | AG-25A-0093 | 100 μg | Rat | ELISA, WB | Ms |
| | ^ | | | | 0 |
| PROTEINS | PID | SIZE | SOURCE | ENDOTOXIN | SPECIES |
| Programulin (human) (roc) | AG 40A 0069V | 10 110 50 110 | HEK303 Colle | <0.01EU/ug | ш., |

| PROTEINS | PID | SIZE | SOURCE | ENDOTOXIN | SPECIES |
|----------------------------|--------------|---------------|--------------|------------|---------|
| Progranulin (human) (rec.) | AG-40A-0068Y | 10 μg 50 μg | HEK293 Cells | <0.01EU/μg | Hu |
| Progranulin (rat) (rec.) | AG-40A-0194 | 10 μg 50 μg | HEK293 Cells | <0.1EU/μg | Rt |

Related Products

| ANTIBODY | PID | SIZE | ISOTYPE/SOURCE | APPLICATION | SPECIES |
|---------------------------------|-------------|--------|----------------|-------------|---------|
| anti-Granulin C (human), pAb | AG-25A-0090 | 100 μg | Rabbit | ELISA, WB | Hu |
| | | | | | |
| PROTEIN | PID | SIZE | SOURCE | ENDOTOXIN | SPECIES |
| Granulin C (human) (rec.) (His) | AG-40A-0129 | 10 μg | E. coli | <1EU/μg | Hu |

The Progranulin Receptor Sortilin

The SORT1 gene is one of the strongest genetic risk factors for Alzheimer's disease and is associated with frontotemporal dementia (FTD). Sortilin controls PGRN trafficking and lysosomal degradation, but PGRN exerts its multiple functions independent of sortilin. Sortilin down-regulation via blocking antibodies, such as Latozinemab Biosimilar, is a key mechanism in increasing PGRN levels suggesting that sortilin is a potential target to correct PGRN reduction, such as that in patients with frontotemporal dementia caused by GRN mutations.

50 µg




Sortilin (human) (rec.) (His)

AG-40B-0229

Binds to human Progranulin (untagged) (Prod. No. AG-40A-0188Y).

FIGURE: Progranulin (human) (rec.) (untagged) (Prod. No. AG-40A-0188Y) is coated on an ELISA plate at 1 μ g/ml overnight at 4°C. Sortilin (human) (rec.) (His) (Prod. No. AG-40B-0229) or a control His-tagged protein (Prod. No. AG-40B-0177) is added (starting at a concentration of 100nM with a twofold serial dilution) during one hour at RT and the interaction is then detected using an anti-His antibody coupled to HRP.

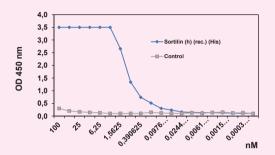


anti-Sortilin (human), mAb (rec.) (blocking) (preservative free) [Latozinemab Biosimilar]

AG-27B-7000PF 100 μg | 500 μg

Clone: Latozinemab **Isotype:** Human IgG1κ

Application: Inhibits sortilin binding to progranulin.

FIGURE: Progranulin (human) (rec.) (untagged) is coated on an ELISA plate at 1 μg/ml. Sortilin (human), mAb (rec.) (blocking) [Latozinemab Biosimilar] or an unrelated mAb (Control) is added (starting at a concentration of 20 μg/ml with a twofold serial dilution) together with 250 ng/ml of Sortilin (rec.) (His) for 1 hour. The decreasing ODs observed in the Y axis represent the binding between Progranulin and Sortilin that is inhibited in a dose-dependent manner by Sortilin (human), mAb (rec.) (blocking) [Latozinemab Biosimilar]. This inhibition is not observed with the control antibody.

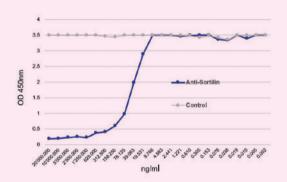

Unique Tag-free Progranulins

Progranulin (human) (rec.) (untagged)AG-40A-0188Y 10 μg | 50 μg

Progranulin (mouse) (rec.) (untagged)AG-40A-0189Y 10 μg | 50 μg

Progranulin (rat) (rec.) (untagged)

AG-40A-0196Y 10 μg | 50 μg

 Higher activity compared to tagged Progranulins (binding to Sortilin)

- Suitable for in vitro and in vivo studies
- Reflects the native sequence with no additional amino acids
- Affinity purified
- Low endotoxin levels (<0.01EU/μg)



TEL +41-61-926-60-40 FAX +41-61-926-60-49 info@adipogen.com

NORTH & SOUTH AMERICA

TEL +1-858-457-8383 FAX +1-858-457-8484 info-us@adipogen.com

Adipogen Corp.

BULK